


1



attach the brush head

- remove brush head from the packet (if not already attached)
- attach by sliding head onto the metal rod of the SmartBrush handle
- store the additional brush head in the packet; replace in 6 weeks

2

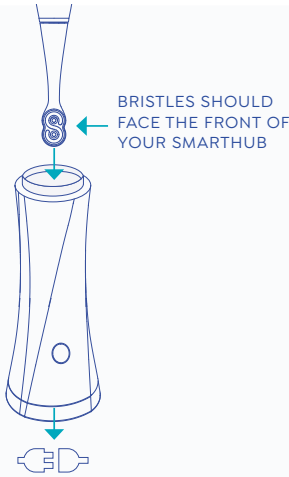

fully charge your devices


your SmartBrush and SmartHub were fully charged at the factory, but may have lost some power in transit

1. remove the cap from the SmartHub
2. place the SmartBrush upside down and insert into the SmartHub; make sure the bristles face the front of your SmartHub
3. plug the charger cable into the wall charger
4. insert the other end into the back of the SmartHub
5. replace the cap, aligning the nodules until it snaps in place

the ring at the base of your SmartBrush will pulse red while charging and turn solid red when fully charged

red light was chosen to avoid disrupting melatonin during nighttime, ensuring a restful sleep

3

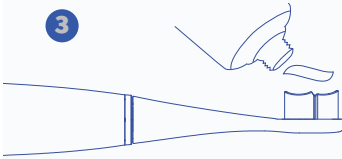


apply toothpaste

- always rinse the brush head
- apply a pea-size amount of toothpaste to the bristles
- put the brush head into your mouth before turning on the brush (to avoid spraying paste around your bathroom!)

for best results, use with SuperMouth MouthPaste

4

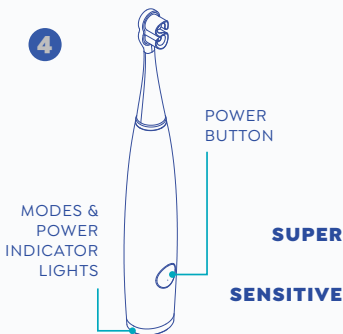


select your preferred mode

- to turn on/off your SmartBrush, press the power button
- the ULTIM8 Sonic SmartBrush has two modes designed to effectively clean your teeth
- to switch modes, turn the toothbrush on and press the button within 10 seconds; the sequence for mode switching is: first mode > second mode > off > first mode

SUPER operates at full power to deliver a vigorous and thorough cleaning. in this mode, the bottom of the handle turns white

SENSITIVE reduced sonic vibrations create a gentler cleaning experience for those with sensitive teeth and gums. in this mode, the bottom of the handle turns blue

5

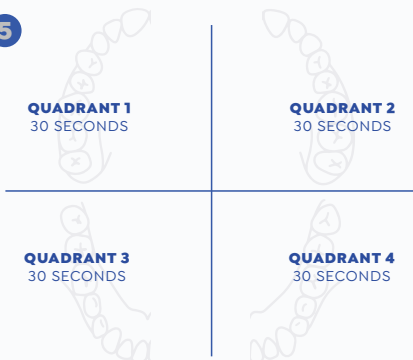

start brushing using the smart timer & pressure sensor

start with the outside surfaces of the teeth. guide the brush head slowly from tooth to tooth, holding the brush head in place for a few seconds against each tooth before moving on to the next one. following along with the shape of each tooth and the curve of the gums is easy with the brush head designed to hug every curve.

about the smart timer

every 30 seconds while brushing, you'll feel a short pause of the vibration, signaling it's time to switch to a different section of your mouth. this ensures each of the four areas receives equal attention. after two minutes, the toothbrush automatically turns off.